



RUBBER, PLASTICS AND ELASTIC MATERIALS DENSITY TESTER TTD-SLPY

Rubber, Plastics and Elastic Materials Density Tester TTD-SLPY is an advanced instrument designed for precise measurement of density and porosity in rubber, plastics, and elastic materials. It is available with maximum weight capacities of 2000g, 3000g, and 6000g, offering a weighing precision of 0.01g (0.1g for the 6000g model) and a density precision of 0.001g/cm³.

Density Meter, Density Tester, Densimeter, Specific Gravity Tester

Density Meter, Density Tester, Densimeter, Specific Gravity Tester

Maximum Weight Capacity: 2000g / 3000g / 6000g

Weighing Precision: 0.01g / 0.01g / 0.1g

Density Precision: 0.001g/cm³

- According to ASTM C20, C128, C127, C134, C135, C437, C357, AASHTOT 84 standard.
- The machine is equipped with 4 functions.
- Show Bulk density, Wet density, Porosity, Absorption, Apparent DS, Open pores, Close porosity, Total porosity directly.

APPLICATIONS

Suitable for

Shoe materials, Sponge, EVA materials, EPDM plastic industry, Graphite, Carbon electrode, Refractory, Corrosion-resistant materials, New materials research lab, Building industry, Civil Engineering industry and concrete mixture ratio research lab, Mineral mining industry, Geology research lab.

Principle

According to ASTM C20, C128, C127, C134, C135, C437, C357, AASHTOT 84 standard and by adopting the buoyancy method of Archimedean principle. and by applying the Vacuum saturation method, the Boiling saturation method can show the measuring result directly.

Mode

1. **BULK** : Adopting the water saturation method for a porous material.
2. **BK+** : Adopting water saturation method for porous absorbent material. With five groups of memory functions.
3. **Powder** : For powder products, use a pycnometer (optional) to directly read the true density of the powder.
4. **DS** : Nonabsorbent material.

FEATURES


- The infrared temperature sensor automatically detects water temperature and performs compensation accordingly.
- Rubber, Plastics and Elastic Materials Density Tester TTD-SLPY equipped with automatic, manual, and solution temperature compensation functions.
- Provides direct readings of body bulk density, wet density, porosity, absorption, apparent density, open pores, closed porosity, and total porosity.
- Includes five groups of memory functions, allowing simultaneous measurement of five sample groups in absorbent product density measurement mode.
- Directly displays values for the development research of new materials. With HI and LO limit functions and a buzzer, the device determines whether the specific gravity of the test object is qualified in non-absorbent product density measurement mode. It also features a function for two groups mixing ratio.
- Offers the option to purchase a pycnometer for direct measurement of powder true density.
- Equipped with RS-232 & USB interface.

- Utilizes a large-size water tank design to minimize errors caused by the buoyancy of the suspension rail.
- The water tank size is 290 × 255 × 160mm.

TECHNICAL SPECIFICATIONS

Model		TTD-SLPY2000	TTD-SLPY3000	TTD-SLPY6000	
Max weight		2000g	3000g	6000g	
Weighting precision		0.01g		0.1g	
Density precision		0.001 g/cm ³			
Wet density		0.001 g/cm ³			
Porosity		0.01 %			
Absorption		0.01 %			
Apparent DS		0.001 g/cm ³			
Open pores		0.001 g/cm ³			
Close porosity		0.01 %			
Total porosity		0.01 %			
Mix ratio :		0.01 %			
Mode	Bulk density (Adopting water saturation method)	Porous material : Show Bulk density, Wet density, Porosity, Absorption, Apparent DS, Open pores, Close porosity, Total porosity			
	Bulk density (Adopting water saturation method) (Five groups of memory functions)	Porous absorbent material : Show Bulk density, Wet density, Porosity, Absorption, Apparent DS, Open pores, Close porosity, Total porosity			
	Powder true density	Powder products : Show Powder true density			

Nonabsorbent material	Show Density, Volume, Mix ratio directly
-----------------------	--

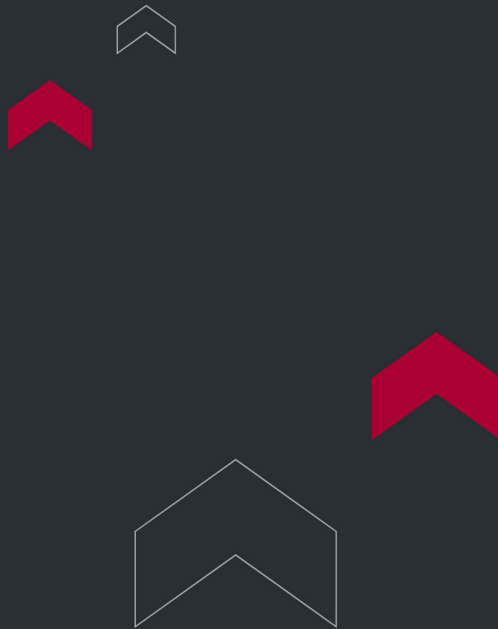

Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 🇨🇦

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 🇺🇸

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 🇦🇪

Phone: +971 488 19252
Email: me@torontech.com

