



FACE MASK DIFFERENTIAL PRESSURE TESTER

Face Mask Differential Pressure Tester ToronFM™ DPT measures the air exchange pressure of medical face masks and other fabrics.

Standard

1. ASTM F2100 Standard Specification for Performance of Materials Used in Medical Face Masks
2. BS EN 14683 Medical face masks - Requirements and test methods, section 5.2.3, Annex C

APPLICATIONS

The Torontech ToronFM™ DPT is a specialized instrument for determining the differential pressure (Delta P) of medical face mask materials, in full compliance with standards like ASTM F2100 and EN 14683. Its primary application is to provide a precise, quantitative measure of a material's breathability. This is a fundamental comfort and performance parameter, directly related to the ease of breathing for the wearer, and is a key criterion for classifying medical masks.

This critical comfort and performance test is a cornerstone of quality and development in:

- **Medical Device Manufacturing:** As an essential quality control and lot release test for manufacturers of medical and surgical face masks, providing the definitive data required to certify their products to the required performance levels (e.g., Type I, II, or IIR) based on breathability.
- **Third-Party & Regulatory Testing:** For independent laboratories and government agencies to verify that face masks being sold in the market meet the breathability requirements of national and international public health and safety standards.
- **Non-Woven & Filter Media R&D:** As a primary research tool for developing new, advanced filter materials, allowing formulators to optimize the balance between high filtration efficiency and low breathing resistance to maximize both protection and user comfort.


- **Hospital & Healthcare Procurement:** To perform incoming quality assurance and supplier validation, ensuring that the masks procured for use by healthcare workers provide the necessary level of comfort and breathability for extended wear during long shifts.

FEATURES

- Powered by a suction air source, offering flexibility in test site locations.
- Equipped with a high-precision pressure difference sensor for digital display of pressure differences.
- Features a specialized sample holder to securely grip the sample for accurate testing.

TECHNICAL SPECIFICATIONS

Model	ToronFM™ DPT
Control mode	7-inch Touch screen
Air source	Suction type
Air flow for testing	8±0.2L / min
Air flow range	0~20L/min
Sealing method	Face seal
Air permeability diameter of sample	Φ25mm
Range of differential pressure sensor	0 ~ 1000Pa, accuracy: ±1%
Weight	40kg
Voltage	Single Phase AC220V 50Hz (110V is also available)

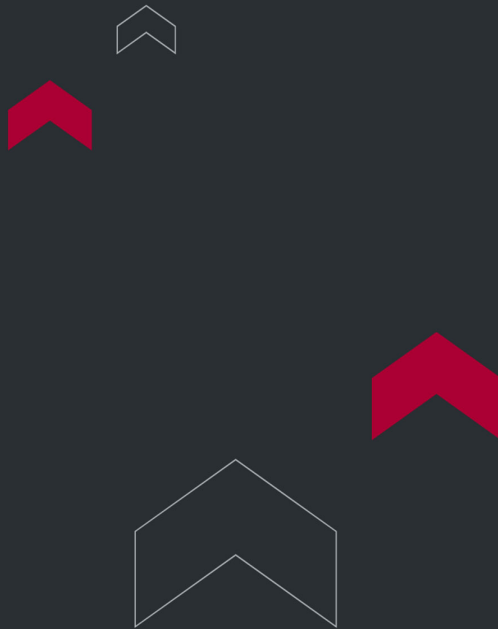

Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 

Phone: +971 488 19252
Email: me@torontech.com

