



TREATMENT CART (PREMIUM STAINLESS)

Treatment Cart (Premium Stainless) - ToronCAB™ 1650

ToronCAB™ 1650 is a wide stainless steel treatment trolley that pairs three open shelves with a removable waste bag for ward rounds. Three railed stainless decks on the left hold supplies and dressings, while a fabric waste bag on the right collects soiled materials during a round. The oxford cloth bag lifts out for washing, which keeps the cart hygienic between rounds. Four 100 mm casters, two with brakes, and push handles at each end move the cart and hold it steady.

Nurses and clinical staff use the Torontech ToronCAB™ 1650 for treatment, dressing, and service rounds that cover several patients. The three shelves carry enough supplies for a full round, and the integrated waste bag collects soiled dressings or linens as staff go, so used items leave the bedside immediately. Full stainless steel construction resists corrosion and tolerates strong disinfectants, while the removable bag washes separately to keep the cart clean. Push handles at each end let staff steer the wide cart from either side.

The wider L950 mm body adds shelf and waste capacity for busy rounds, yet stays mobile through corridors on its braked casters. Built for daily use, the ToronCAB™ 1650 brings high-capacity supply storage and integrated waste collection together in one durable, fully cleanable cart.

APPLICATIONS

Treatment Cart (Premium Stainless) - ToronCAB™ 1650 Applications

The ToronCAB™ 1650 supports treatment, dressing, and service rounds where staff carry supplies and collect waste across several patients. Its three shelves and removable waste bag fit a range of ward and clinical tasks.

- Treatment and Dressing Rounds: carries supplies on three shelves and collects soiled dressings in the waste bag while staff move through a ward.
- Hospital Inpatient Wards: serves as a mobile supply-and-waste station for routine care across multiple beds.

- Housekeeping and Linen Collection: gathers soiled linens or packaging in the removable bag during cleaning and changeover rounds.
- Outpatient and Day Clinics: restocks treatment bays and removes waste between patients.
- Operating Room and Recovery Support: moves supplies and collects waste during turnover between cases.
- Long-Term Care and Rehabilitation Facilities: supports routine rounds with supply storage and waste collection in one cart.
- Central Sterile and Supply Departments: transports supplies and collects packaging or waste during distribution.



FEATURES

Treatment Cart (Premium Stainless) - ToronCAB™ 1650 Key Features

The ToronCAB™ 1650 combines three-tier supply storage with integrated, removable waste collection in a wide, all-stainless treatment cart.

- Three Railed Stainless Shelves: three open decks with siderails hold a full round of supplies, instruments, and dressings.
- Removable Fabric Waste Bag: a large oxford cloth bag collects soiled dressings, linens, or packaging during a round and lifts out for washing.
- Full Stainless Steel Construction: a corrosion-resistant frame tolerates strong disinfection in clinical and sterile areas.
- Dual Push Handles: handles at each end let staff steer the wide cart from either side.
- Four 100 mm Casters with Brakes: four casters roll smoothly over hospital floors, and two with brakes hold the cart steady when parked.
- Wide, High-Capacity Build: the L950 mm body adds shelf and waste capacity for busy multi-patient rounds.

THEORY & METHOD

Theory and Method


The ToronCAB™ 1650 splits a round into two jobs on one frame: carrying clean supplies and removing waste. Three open shelves on one side give staff visible storage for dressings, instruments, and consumables, with a siderail on each to stop items sliding off. A large fabric waste bag on the other side takes soiled materials as the round progresses, so used items leave each bedside without a separate trip to a bin. Keeping supply and waste together lets one cart handle a full round.

Construction keeps the cart clean and easy to move. Full stainless steel resists corrosion and tolerates strong disinfection, while the waste bag lifts out of its frame to wash separately rather than soiling a fixed bin. Push handles at each end let staff steer the wide cart from either side or guide it through a doorway. Four 100 mm casters carry the loaded cart over hospital floors, with brakes on two to hold it steady while staff load supplies or empty the bag.

TECHNICAL SPECIFICATIONS

Treatment Cart (Premium Stainless) - ToronCAB™ 1650 Technical Specifications

Parameter	Specification
External Dimensions	950 (L) × 500 (W) × 900 (H) mm
Construction	Full stainless steel
Shelves	Three stainless steel decks with siderails (left side)
Waste Collection	Removable oxford cloth waste bag (right side)
Handle	Lateral push handle at each end
Casters	Four × 100 mm; two with brakes
Packing Size	990 × 770 × 560 mm

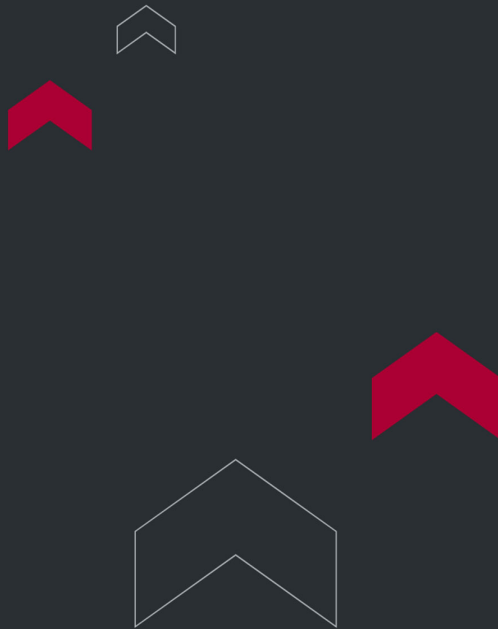

Torontech Inc.

251 Consumers Road, Suite 1200
Toronto, Ontario, M2J 4R3, Canada 

Phone: +1 416 368 2721 | Fax: +1 416 981 7652
Email: info@torontech.com

Torontech, LLC

601 Brickell Key Drive, Suite 700
Miami, FL 33131, USA 

Toll-Free: 1-866-383-7919
Email: sales@torontech.com

Torontech ME FZE

Meydan Grandstand, 6th Floor,
Meydan Road, Nad Al Sheba,
Dubai-261440, UAE 

Phone: +971 488 19252
Email: me@torontech.com

